

ORION®

INSTRUCTION MANUAL

INFRARED MASSAGER
Model: OR-MS23

ENGLISH

GB



Dear Customer! Thank you for buying ORION products. Carefully read all instructions before operating and save this instruction manual for future reference.

<http://www.orion.ua>

**Dear Customer, Thank you for buying ORION products!
Welcome to use!**

IMPORTANT SAFEGUARDS

Before connecting the appliance check that the voltage indicated on the type plate corresponds to the mains voltage in your country.

When using electrical appliances basic safety precautions should always be followed. Please read the following carefully before use:


1. Always unplug the appliance immediately after use.
2. Do not use outdoor.
3. Do not use while in the bath.
4. Do not use with wet hands.
5. Do not allow the appliance to be placed or dropped into water or any other liquid.
6. Individuals with pacemakers should not use the appliance.
7. Do not use after taking pain killers, sedatives or alcohol; your sensitivity to heat could be reduced.
8. The appliance is designed to provide only temporary relief for minor aches and pains muscular tensions or similar and improve your well-being. Do not use the appliance as a substitute for medical attention.
9. Do not reach for an appliance that has fallen into water. Turn off the power supply to the appliance and then unplug from the wall socket.
10. Do not leave the appliance unattended when it is switched on.
11. Always allow the appliance to cool down before packing away.
12. Close supervision is necessary when this appliance is used by, on, or near children or invalids.
13. The appliance should be used for its intended purpose only.
14. Never use attachments not recommended or sold by the appliance manufacturer.
15. Keep the cord away from heated surfaces.
16. Check the cord regularly for wear or damage, particularly where the cord enters the appliance and the plug.
17. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Call the Consumer Service Centre for examination and repair.
18. Do not wrap the power cord around the appliance, straighten after use.
19. Be careful when placing hot surfaces on areas sensitive to heat such as the face and neck.
20. Do not use on or near the eyes.
21. Take special care when using on children or anyone who is sensitive to heat.
22. Switch the appliance off before putting it down, even if it is just for "a moment".
23. Do not allow the air inlet to become blocked or even partially blocked. This may cause the appliance to overheat.
24. Do not place any objects into the air inlets or any other openings.
25. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliances by a person responsible for their safety.
26. Children should be supervised to ensure that they do not play with the appliance.

CAUTION!

Please take note of the following when using the Infrared and Magnetic settings:

1. Do not use for more than 20 minutes on a specific area.
2. Do not use if you have sensitive skin, poor circulation, swelling, inflammation, or sores are present.
3. Do not use if you are experiencing unexplained pain or swollen muscles, or following a serious muscular injury before consulting your physician.
4. Discontinue use of massager if pain increases or swelling occurs after application.
5. Do not use in the presence of phlebitis and thrombosis, open and fresh wounds, bruises, broken skin or varicose veins.
6. Do not use the vibrating massager directly on one concentrated area of the body. Vibrating massagers are designed to be used in a continuous motion, constantly changing the massage area.
7. If you have diabetes, poor circulation or a heart condition, consult your physician prior to use.

PARTS

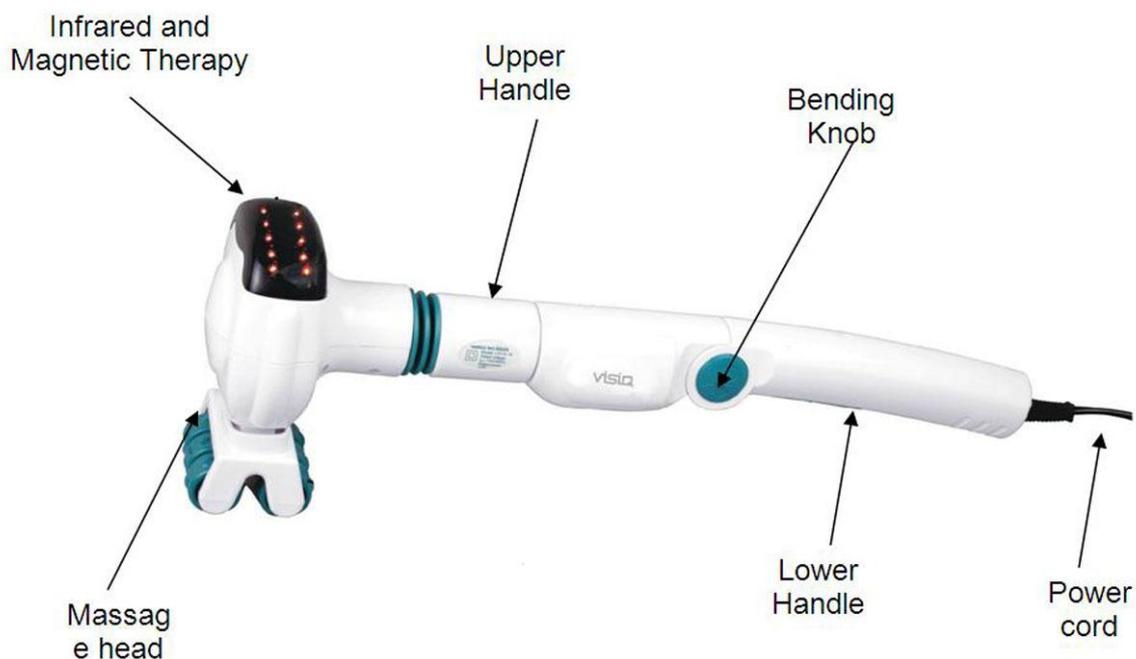

To increase the effectiveness of massage we recommend to use special attachments. Four interchangeable attachments provide a full body massage: back, shoulders, arms, legs, hips and neck.

You can massage the body through the clothing, but the best effect is achieved while massaging a naked body.

It is best to attach the desired massage attachment when the massager is off.

For massaging sensitive areas of the body, as well as using infrared heating we recommend to use the device without attachments.

Massage Attachments:

1. Rolling Massage Head
2. Body Massage Head
3. Joint Massage Head
4. Concentrator Head



HOW TO USE

Use of this Infrared massager has a relaxing, soothing and stimulating effect on your muscles. If used regularly for several minutes a day, it can reduce stress, fatigue and muscular pain, and it has anti-cellulite function.


This appliance is equipped with an infrared lamp. The infrared radiation penetrates into the deep layers of skin causing the pores and vessels to expand and simultaneously enhancing the metabolism process and microcirculation. Using the infrared massager brings relief from muscular pain and regenerates the skin.

How to use:

Hold the handle of the device and select the desired function in accordance with the marks on the appliance.

Move the massager in the directions indicated on the picture.

Note: these directions are recommended rather than compulsory.

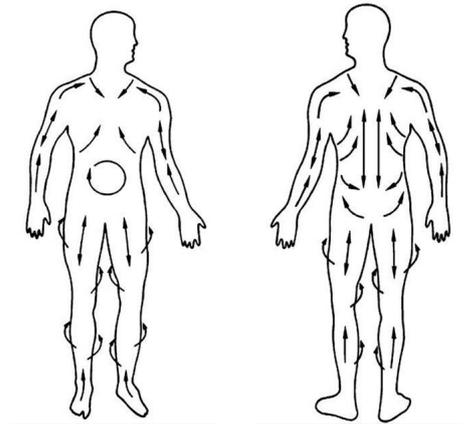

OFF (0)	The appliance is off.
Massage (M)	This provides massage only. Best for use with your massage attachments.
MH	This setting switches on the Infrared and Magnetic Therapy, and provides massage. This setting is most effective without the use of attachments.

CARE AND MAINTENANCE

- Switch the appliance off, unplug it and let it cool down.
- Wipe the body of the massager with a soft damp cloth. If necessary use a little washing-up liquid.
- Remove deposits from the inlet grille by brushing with a soft brush.
- Do not immerse the appliance in water or any other liquid.
- Do not use harsh or abrasive cleaning agents or solvents.

TECHNICAL CHARACTERISTICS

- Voltage: 230V
- Power: 12W
- Rated current: 0,05A
- Water proof: IPX0
- Net/gross weight: 690g/1050g
- Product size: 420X105X95mm
- Gift box size: 480X105X115mm
- Colour: white/green